Rat Intact PTH ELISA Kit

Enzyme-Linked ImmunoSorbent Assay (ELISA) for the Quantitative Determination of Rat Intact Parathyroid Hormone Levels in Serum, Plasma or Cell Culture Media

Immutopics

Immutopics, Inc.

96 Test Kit Cat.# 60-2500

Store at 2 - 8°C Upon Receipt

For RESEARCH Use Only

Not for use in diagnostic procedures

INTENDED USE

This kit is intended for research use only in the determination of rat intact PTH levels in serum, plasma or cell culture media.

INTRODUCTION

Rat intact parathyroid hormone (PTH) is an 84 amino acid polypeptide produced by the parathyroid gland with its biological activity residing in the N-terminal region of the peptide. PTH plays an important role in maintaining the concentration of ionized calcium within the limits needed to achieve normal metabolic functions. When serum calcium levels are decreased the parathyroid gland increases secretion of the hormone which results in increased mobilization of calcium from skeletal reserves into the circulation. When levels of serum calcium are increased the secretion of PTH is reduced.

The similarities between rat and human physiology relative to calcium metabolism make the rat an excellent live-animal model for studying skeletal disease and in the pre-clinical evaluation of pharmacologic agents that may alter bone remodeling. Quantitation of the biologically active rat intact PTH with this kit can provide a precise and sensitive assessment of changes in bone and mineral metabolism.

TEST PRINCIPLE

The Rat Intact PTH ELISA Kit is a two-site enzyme-linked immunosorbent assay (ELISA) for the measurement of PTH in rat serum, plasma or cell culture media. Two different goat polyclonal antibodies to rat intact PTH have been purified by affinity chromatography. The antibody which recognizes epitopes within the midregion/C-terminal portion (39-84) of the peptide is biotinylated for capture. The other antibody which recognizes epitopes within the N-terminal region (1-34) is conjugated with the enzyme horseradish peroxidase (HRP) for detection

A sample containing rat intact PTH is incubated simultaneously with the biotinylated capture antibody and the HRP conjugated antibody in a streptavidin coated microtiter well. Intact PTH (1-84) contained in the sample is immunologically bound by the capture antibody and the detection antibody to form a "sandwich" complex:

Well/Avidin-Biotin Anti-Rat PTH — Rat Intact PTH — HRP Anti-Rat PTH

At the end of this incubation period, the well is washed to remove any unbound antibody and other components. The enzyme bound to the well is then incubated with a substrate solution in a timed reaction and then measured in a spectrophotometric microtiter plate reader. The enzymatic activity of the antibody complex bound to the well is directly proportional to the amount of rat intact PTH in the sample. A standard curve is generated by plotting the absorbance versus the respective rat intact PTH concentration for each standard on linear or logarithmic scales. The concentration of rat intact PTH in the samples is determined directly from this curve.

REAGENTS: Preparation and Storage

Store the kit at $2-8^{\circ}$ C upon receipt. Store the standards and controls at -20° C or below after reconstitution. For the expiration date of the kit refer to the label on the kit box. All components are stable until this expiration date.

Prior to use allow all reagents to come to room temperature and mix by gentle swirling and inversion. Reagents from different kit lot numbers should not be combined or interchanged.

1. STREPTAVIDIN COATED MICROTITER PLATE (40-0010)

One plate with 12 eight well strips and frame (96 wells total). This reagent should be stored in the foil pouch with desiccant at 2 - 8°C and is stable until the expiration date on the kit.

2. RAT PTH BIOTINYLATED ANTIBODY (40-2510)

One vial containing 5.5 mL of biotin labeled anti-rat PTH in TRIS buffered saline with protein stabilizers and a non-azide, non-mercury preservative. This reagent should be stored at 2 - 8°C and is stable until the expiration date on the kit.

3. RAT PTH HRP CONJUGATED ANTIBODY (40-2520)

One vial containing 5.5 mL of horseradish peroxidase conjugated to anti-rat PTH in a stabilized protein solution with a non-azide, non-mercury preservative. This reagent should be stored at 2 - 8°C protected from light and is stable until the expiration date on the kit.

NOTE: Make a Working Antibody Solution by pipetting equal volumes of Rat PTH Biotinylated Antibody and Rat PTH HRP Conjugated Antibody prior to use. Mix only the volume required for immediate use. Mix well to ensure homogeneity.

4. RAT INTACT PTH STANDARDS (40-2531 to 40-2536)

Six vials each containing rat intact PTH (1-84) lyophilized in a protein matrix with a non-azide, non-mercury preservative. **Refer to vial label for exact concentration.** Before use reconstitute the vial with the rat intact PTH concentration of 0 pg/mL with 2.0 mL of deionized water. Before use reconstitute each of the other five vials of standards with 1.0 mL of deionized water. Allow the vials to sit for approximately 20 minutes with occasional gentle swirling and inversion. Assure complete reconstitution before use.

Use the standards immediately after reconstitution; freeze the unused portion for later use. After reconstitution the standards are stable until the expiration date on the kit box when stored at -20°C or below with up to 3 freeze/thaw cycles.

5. RAT INTACT PTH CONTROLS I & II (40-2541 & 40-2542)

Two vials each containing rat intact PTH (1-84) lyophilized in a protein matrix with a non-azide, non-mercury preservative. **Refer to vial label for control ranges.** Before use reconstitute each control with 1.0 mL of deionized water. Allow the vials to sit for approximately 20 minutes with occasional gentle swirling and inversion. Assure complete reconstitution before use.

Use the controls immediately after reconstitution; freeze the unused portion for later use. After reconstitution the controls are stable until the expiration date on the kit box when stored at -20°C or below with up to 3 freeze/thaw cycles.

6. ELISA WASH CONCENTRATE (40-0041)

One vial containing 20 mL of a 20 fold concentrate. Before use dilute the contents to 400 mL with deionized water and mix well. Upon dilution this yields a working wash solution containing a surfactant in phosphate buffered saline with a non-azide, non-mercury preservative. The diluted wash solution should be stored at room temperature and is stable until the expiration date on the kit.

7. ELISA HRP SUBSTRATE (40-0027)

One bottle containing 16 mL of tetramethylbenzidine (TMB) with hydrogen peroxide. This reagent should be stored at 2 - 8°C protected from light and is stable until the expiration date on the kit.

8. ELISA STOP SOLUTION (40-0030)

One bottle containing 11 mL of 1 M sulfuric acid. This reagent may be stored at room temperature or at 2 - 8° C and is stable until the expiration date on the kit.

9. PLATE SEALER (10-2016)

Two included in kit.

SAFETY PRECAUTIONS

Avoid contact with reagents containing TMB, hydrogen peroxide, or sulfuric acid (i.e. ELISA HRP Substrate and ELISA Stop Solution). In case of contact with any of these reagents, wash thoroughly with water. TMB is a suspected carcinogen. Use Good Laboratory Practices. Wash hands before eating. Do not eat, drink or smoke in the work area.

MATERIALS REQUIRED BUT NOT PROVIDED

- 1.0 mL and 2.0 mL volumetric pipettes for reconstituting standards and controls.
- 2. Precision pipets capable of delivering 25 μL ,100 μL and 150 μL.
- 3. Aluminum foil.
- 4. Automated microtiter plate washer OR
- Repeating dispenser for delivering 350 μL and suitable aspiration device.
- Container for storage of wash solution.
- Spectrophotometric microtiter plate reader capable of reading absorbance at 450 nm and at 595 - 650 nm.
- 8. Deionized water.
- 9. Horizontal rotator capable of maintaining 180 220 RPM.
- 10. Timer.

SPECIMEN COLLECTION

The measurement of the rat intact PTH concentration may be made using serum, plasma or cell culture media. Fifty microliters of serum, plasma, or culture media are required to assay the sample in duplicate. If obtaining serum, collect blood and allow it to clot at room temperature. Centrifuge the sample and separate the serum, plasma or media from the cells. Samples should be assayed immediately or stored frozen at -20°C or below. Avoid repeated freezing and thawing of specimens.

The use of various anesthetics can cause significant elevations in serum and plasma PTH concentrations. It is therefore imperative to use consistent sample collection procedures within studies. (See ref. #4 and #6)

ASSAY PROCEDURE

- Place a sufficient number of Streptavidin Coated Strips in a holder to run PTH standards, controls and unknown samples.
- Pipet 25 μL of standard, control, or sample into the designated or mapped well. Freeze the remaining standards and controls as soon as possible after use.
- Pipet 100 μL of the Working Antibody Solution consisting of equal volumes of Rat PTH Biotinylated Antibody and Rat PTH HRP Conjugated Antibody into each well.
- Cover the plate with one plate sealer, then cover with aluminum foil to avoid exposure to light.
- Incubate plate at room temperature for three hours on a horizontal rotator set at 180 - 220 RPM.
- 6. Remove the aluminum foil and plate sealer. Using an automated microtiter plate washer aspirate the contents of each well. Wash each well five times by dispensing 350 μL of working wash solution into each well and then completely aspirating the contents. A suitable aspiration device may also be used.

- Pipet 150 μL of ELISA HRP Substrate into each of the wells.
- Re-cover the plate with the Plate Sealer and aluminum foil. Incubate at room temperature for 30 minutes on a horizontal rotator set at 180 - 220 RPM.
- Remove the aluminum foil and plate sealer. Read the absorbance at 620 nm (see Note) within 5 minutes in a microtiter plate reader against the 0 pg/mL Standard wells as a blank.
- Immediately pipet 100 μL of ELISA Stop Solution into each of the wells. Mix on horizontal rotator for 1 minute.
- 11. Read the absorbance at 450 nm within 10 minutes in the microtiter plate reader against a reagent blank of 150 μ L of Substrate and 100 μ L of Stop Solution.

If dual wavelength correction is available set the Measurement wavelength to 450 nm and Reference wavelength to that used in step #9.

NOTE: Absorbance may be read at wavelengths from 595 nm to 650 nm depending upon available filters.

PROCEDURAL NOTES

- It is recommended that all standards, controls and samples be assayed in duplicate. The average absorbance reading of each duplicate should then be used for data reduction and the calculation of results.
- Store light sensitive reagents (i.e. HRP Conjugated Antibody, the Working Antibody Solution consisting of combined Biotinylated Antibody and HRP Conjugated Antibody, and ELISA HRP Substrate) in the original amber bottles or other suitable container which is well protected from light.
- 3. Store any unused Streptavidin Coated Strips in the resealable aluminum pouch with desiccant to protect them from moisture.
- The sample and all reagents should be pipetted carefully to minimize air bubbles in the wells.
- 5. The sequence and timing of each reagent addition is important as both the immunological and enzymatic reactions are in kinetic modes. The washing step is also an important part of the total assay procedure. The use of an automated microtiter plate washer is strongly recommended. All pipeting and washing steps should be performed such that the timing is as consistent as possible.
- Samples with values greater than the highest standard should be diluted 1:10 with the 0 pg/mL Standard or a suitable sample diluent and reassayed. Multiply the result by 10. (See Limitations, # 1 and # 2)
- 7. Plasma or cell culture media samples may contain fibrin clots or cellular debris. Freeze/thaw of plasma samples may accelerate clot formation. These samples must be centrifuged and decanted prior to assay to remove all particulate material which can cause random high non-specific binding on well surface.
- Rarely, upon opening the streptavidin plate, small white crystals
 may be observed in some of the wells. This is entirely cosmetic
 and will not affect the assay. This condition is reported by other
 kit manufacturers and results from the final stabilizing buffer used
 in the coating process.

CALCULATION OF RESULTS

The two absorbance readings taken before and after the addition of the ELISA Stop Solution allow for the construction of two standard curves using the rat intact PTH standards contained in the kit. **Refer to the individual vial label for exact concentration.** The primary curve used for calculation of results is the second reading taken after the addition of the ELISA Stop Solution and read at 450 nm. This data utilizes the absorbance values obtained with the first five standards. The first reading taken before the addition of the ELISA Stop Solution and read at 595 nm – 650 nm is intended to extend the analytical range to the value of the sixth (highest) standard provided in the kit. **It should be utilized only for sample results that fall between the value of the fifth and sixth standard.** Results obtained with the first reading should not replace the on-scale reading at 450 nm. Each curve should be generated as follows:

Primary Procedure — Read at 450 nm

- Calculate the average absorbance for each pair of duplicate assay wells.
- Subtract the average absorbance of the 0 pg/mL Standard from the average absorbance of all other readings to obtain corrected absorbance.
- The standard curve is generated by plotting the corrected absorbance of the first five standard levels on the ordinate against the standard concentration on the abscissa using linearlinear or log-log paper. Appropriate computer assisted data reduction programs may also be used for the calculation of rat intact PTH results.

The rat intact PTH concentration of the controls and samples are read directly from the standard curve using their respective corrected absorbance. If log-log graph paper or computer assisted data reduction programs utilizing logarithmic transformation are used, samples having corrected absorbance between the 0 pg/mL Standard and the next highest standard should be calculated by the formula:

Value of unknown = $\frac{\text{(unknown)}}{\text{(unknown)}} x \text{ Value of the 2}^{\text{nd}} \text{ Std.}$ $\frac{\text{Corrected Absorbance}}{\text{(2}^{\text{nd}} \text{ Std.)}}$

Secondary Procedure — Read at 595 nm - 650 nm


- Calculate the average absorbance for each pair of duplicate assay wells.
- The standard curve is generated by plotting the absorbance of the three highest standards on the ordinate against the standard concentration on the abscissa using linear-linear or log-log graph paper.
- The rat intact PTH concentration of samples reading only between the fifth and sixth standard are read directly from this standard curve.

EXAMPLE DATA AND STANDARD CURVE

The following are representative examples of data and the resulting standard curves from the primary and secondary procedures. These curves should not be used in lieu of a standard curve run with each assay.


WELL		ARY ASSAY - AVERAGE	- 450 nm CORRECTED	RESULTS	
I.D.	ABS	ABS	ABS	pg/mL	_
Reagent Blank	0.000 0.000	0.000			
0 pg/mL	0.009 0.011	0.010	0.000		
25 pg/mL	0.062 0.066	0.064	0.054		
80 pg/mL	0.186 0.178	0.182	0.172		
277 pg/mL	0.615 0.603	0.609	0.599		
840 pg/mL	1.841 1.813	1.827	1.817		
Control I	0.117 0.131	0.124	0.114	53	
Control II	0.366 0.366	0.366	0.356	164	
Sample 1	0.516 0.520	0.518	0.508	234	
Sample 2	0.048 0.045	0.047	0.037	17	
Sample 3	2.204 2.200	2.202	2.192	*	
* > 840 pg/mL.	Calculate	using second	ary assay.		

Rat Intact PTH ELISA Primary Assay

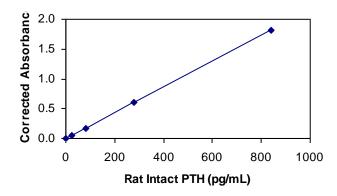

WELL I.D.	SECONDARY ABS	ASSAY - 620 ni AVERAGE ABS	m RESULTS pg/mL
0 pg/mL	0.000	0.000	
277 pg/mL	0.136 0.128	0.132	
840 pg/mL	0.411 0.409	0.410	
2700 pg/mL	1.309 1.322	1.315	
Sample 3	0.491 0.495	0.493	1,012

Rat Intact PTH ELISA Secondary Assay

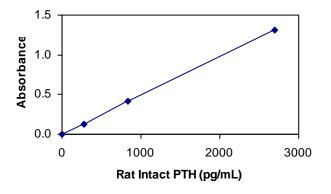

LIMITATIONS OF THE PROCEDURE

- The lowest concentration of rat intact PTH measurable is 1.6 pg/mL (assay sensitivity) and the highest concentration of rat intact PTH measurable without dilution is the value of the highest standard.
- 2. The reagents in this Rat Intact PTH ELISA kit have been optimized so that the high dose "hook effect" is not a problem for samples with elevated rat intact PTH values. Samples with rat intact PTH levels between the highest standard and 200,000 pg/mL will read greater than the highest standard and should be diluted 1:10 with the 0 pg/mL Standard or a suitable sample diluent and reassayed for correct values.
- Grossly lipemic serum or plasma samples may affect the immunological response and it is recommended that results obtained with such samples be scrutinized accordingly.

4. Differences in protein concentration and protein type between samples and standards in an immunoassay contribute to "protein effects" and dose biases. When measuring low protein concentration culture media samples against high protein concentration standards, it is recommended that like samples be assayed together in the same assay to minimize this bias.

QUALITY CONTROL

To assure the validity of the results each assay should include adequate controls with known levels of rat intact PTH. Immutopics recommends that all assays include the laboratory's own rat intact PTH controls in addition to those provided with this kit.

PERFORMANCE CHARACTERISTICS: SENSITIVITY

The sensitivity of the rat intact PTH assay as determined by the 95% confidence limit on 20 duplicate determinations of the 0 pg/mL Standard is 1.6 pg/mL.

PRECISION

To assess intra-assay precision the mean and coefficient of variation were calculated from 20 duplicate determinations of two samples each performed in a single assay.

Mean Value (pg/mL)	Coefficient of Variation
54	2.4 %
173	2.1 %

To assess inter-assay precision the mean and coefficient of variation were calculated from duplicate determinations of two samples performed in 20 assays.

Mean Value (pg/mL)	Coefficient of Variation		
53	6.0 %		
164	5.1 %		

PARALLELISM

Rat serum samples were diluted with the 0 pg/mL Standard and assayed. Results in pg/mL are as follows:

SAMPL	E DILUTION	OBSERVED VALUE	EXPECTED VALUE	% O/E
1	undiluted	124		
	1:2	62	62	100
	1:4	28	31	90
	1:8	15	16	94
2	undiluted	480		
	1:2	233	240	97
	1:4	108	120	90
	1:8	52	60	87
3	undiluted	945		
	1:2	488	473	103
	1:4	231	236	98
	1:8	106	118	90

RECOVERY

Various amounts of rat intact PTH were added to three different rat serum samples and assayed. Results in pg/mL are as follows:

SAMPLE	ORIG. VALUE	AMOUNT ADDED	OBSERVED VALUE	EXPECTED VALUE	% O/E
1	51	212 425 638	261 444 676	263 476 689	99 93 98
2	8	230 460 689	250 431 727	238 468 697	105 92 104
3	0	229 457 686	237 428 790	229 457 686	104 94 115

CROSS-REACTIVITY

The rat intact PTH ELISA is specific for the rat PTH 1-84 molecule. N-terminal 1-34 or mid and C-terminal 39-84 fragments will not be measured. Cross-reactivity with human PTH 1-84 is 2.3%.

WARRANTY

This product is warranted to perform as described in its labeling and literature when used in accordance with all instructions. Immutopics, Inc. DISCLAIMS ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, and in no event shall Immutopics, Inc. be liable for consequential damages. Replacement of the product or refund of the purchase price is the exclusive remedy for the purchaser. This warranty gives you specific legal rights and you may have other rights which vary from state to state.

REFERENCES

- Schultz VL, Garner SC, Lavigne JR, Toverud SU, "Determination of bioactive rat parathyroid hormone (PTH) concentrations in vivo and in vitro by a 2-site homologous immuno-radiometric assay", *Bone* and *Mineral* 27 (1994) 121-132.
- Jara A, Bover J, Lavigne J, Felsenfeld A, "Comparison of Two Parathyroid Hormone Assays for the Rat: The New Immunoradiometric and the Older Competitive Binding Assay", J Bone Miner Res, Vol. 9, No. 10, 1994, p. 1629.
- Rucinski B, Mann GN, Epstein S "A new rapid and reproducible homologous immunoradiometric assay for amino-terminal parathyroid hormone in the rat," Calcif Tissue Int, 1995, 56:83-87.
- Schultz VL, Boass A, Garner SC, Toverud SU, "Several Anesthetics, but Not Diethyl Ether, Cause Marked Elevation of Serum Parathyroid Hormone Concentration in Rats", *J Bone Miner Res*, Vol. 10, No. 9, 1995, p. 1298.
- Urena P, Kubrusly M, Mannstadt M, Hruby M, et al., "The renal PTH/PTHrP receptor is down-regulated in rats with chronic renal failure", Kidney International, Vol. 45 (1994), pp. 605-611.
- Mallya SM, Buchberger G, Arnold A, "Effects of anesthetic agents on serum parathyroid hormone and calcium concentrations in mice." Veterinary Anesthesia and Analgesia, 2007, 34: 403-407.

CLIENT SERVICES

To place an order or for technical assistance, contact Immutopics International at (800) 681-6665 or (949) 369-9207 or FAX to (949) 369-9405 or e-mail: clientservices@immutopicsintl.com.

Developed and Immutopics, Inc.
Manufactured by: San Clemente, CA 92673

Distributed by: Immutopics International
San Clemente, CA 92673

www.immutopicsintl.com

Catalog # 60-2500

90-2500

Effective: 07/14

Copyright© *Immutopics, Inc.*San Clemente, CA