

## Information

*always your partner*



# TECO<sup>®</sup> Fast Cryptococcus neoformans Antigen LFA


For early diagnosis of Cryptococcus infections such as Cryptococcus meningitis, and as a monitoring tool during therapy

- Fluorescence-based lateral flow assay designed for readout via reader
- Quantitative determination of the Cryptococcus antigen glycuronoxylomannan (GXM) in ng/mL in only 20 minutes
- Detects the respective serotypes A, B, C, D and AD
- Clinically validated for CSF and serum samples
- No sample treatment required
- Limit of detection: 0.1 ng/mL, clinical cut-off value: 0.5 ng/mL

	Serum (n = 310)	CSF (n = 350)
Clinical sensitivity	97%	100%
Clinical specificity	98.1%	100%



## TECO® *Fast Cryptococcus neoformans* Antigen LFA Procedure

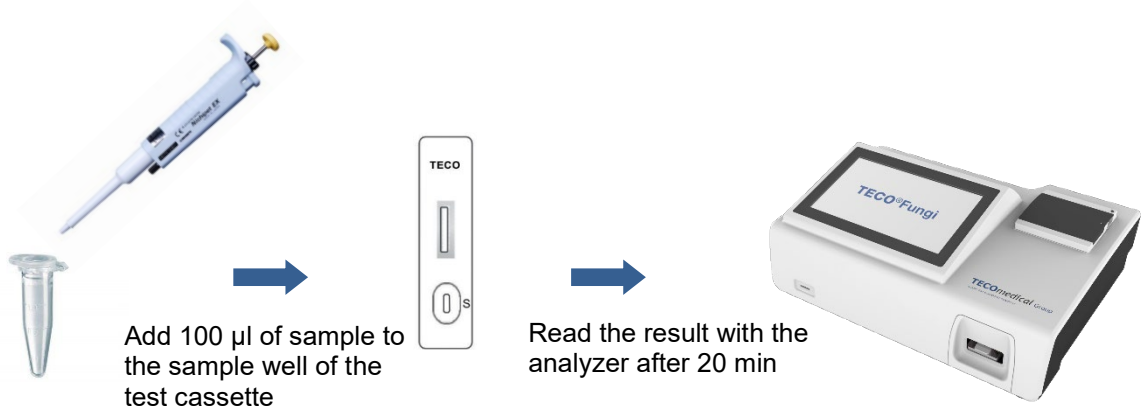

TECO® *Fast Cryptococcus neoformans* Antigen LFA  
Catalogue No.: TE1085

Für weitere Informationen wenden Sie sich bitte an:  
For further information please contact:  
Pour plus d'information, veuillez contacter:

[www.tecomedical.com](http://www.tecomedical.com)



A EUROBIO SCIENTIFIC COMPANY

Switzerland / Headquarters  
**TECOmedical AG**  
Gewerbstrasse 10  
4450 Sissach  
Phone +41 61 985 81 00  
Fax +41 61 985 81 09  
Mail [info@tecomedical.com](mailto:info@tecomedical.com)

Germany  
**TECOmedical GmbH**  
Wasserbreite 57  
32257 Bünde  
Phone +49 52 23 985 99 99  
Fax +49 52 23 985 99 98  
Mail [info@tecomedical.com](mailto:info@tecomedical.com)

Austria  
**TECOmedical AG**  
Phone 0800 20 40 66  
Fax 0800 20 40 55  
Mail [info@tecomedical.com](mailto:info@tecomedical.com)