

Information

always your partner

Real Accurate Quadruplex FLU/COVID-19 RT PCR




RT-PCR for extraction-free determination of Influenza and COVID-19

- Detection of Influenza A, Influenza B and SARS-CoV-2
- Validated, simple treatment protocol for thermal sample preparation
 - Time and cost savings (efficiency)
 - Independent from availability of extraction material
 - Discharge of extraction systems
- CE certified IVD method (RT-PCR plus thermal treatment)

Test principle	Detection of viral RNA from Influenza A, Influenza B and SARS-CoV-2 in clinical samples
Technology	Real-time RT-PCR
Sample type	Nasopharyngeal swab
Controls included	Positive control, negative control and internal control
Validated on	LC480®, Rotor-Gene® Q, Quantstudio 5®, CFX96®, Mic®, MyGoPro ESR®
LoD (µl/reaction/Thermal treatment)	Flu A (H3N2) : 5.5 ; Flu A (H1N1 pdm09) :14.5 ; Flu B : 28.0; SARS-CoV-2 : 5.5.
Specificity	No crossreactivity in samples with 34 relevant pathogens
Thermal treatment	10 minutes at 99°Celsius, Centrifugation, 5µl into PCR mix
Validated Media	E-Swab Media , VTM and 0.9% (w/v) NaCl
Description	RealAccurate Quadruples Flu/COVID-19 PCR Kit
Produced by	PathoFinder B.V., Maastricht, The Netherlands
Packaging	50 reactions
Reference	PF0976-R

Sample preparation using the thermal treatment protocol

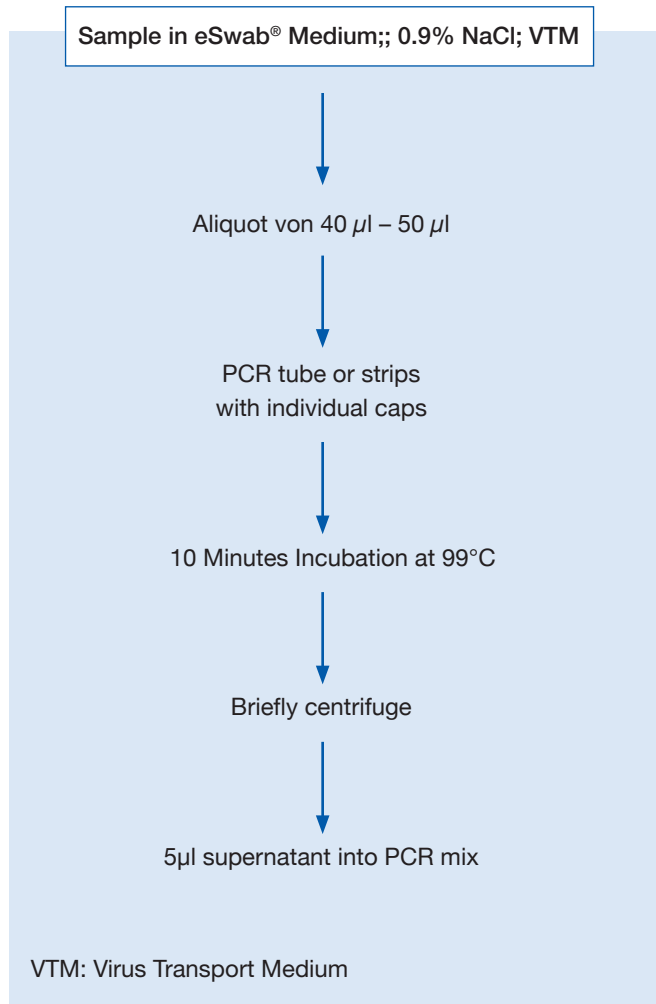

Für weitere Informationen wenden Sie sich bitte an:
For further information please contact:
Pour plus d'information, veuillez contacter:

www.tecomedical.com



A EUROBIO SCIENTIFIC COMPANY

Switzerland / Headquarters

TECOmedical AG

Gewerbestrasse 10

4450 Sissach

Phone +41 61 985 81 00

Fax +41 61 985 81 09

Mail info@tecomedical.com

Germany

TECOmedical GmbH

Wasserbreite 57

32257 Bünde

Phone +49 52 23 985 99 99

Fax +49 52 23 985 99 98

Mail info@tecomedical.com

Benelux

TECOmedical Benelux BV

Computerweg 22

3542 DR Utrecht, The Netherlands

Phone +31 30 307 87 30

Fax +31 30 307 49 39

Mail benelux@tecomedical.com

Austria

TECOmedical AG

Phone 0800 20 40 66

Fax 0800 20 40 55

Mail info@tecomedical.com